# HOWTO LED

(Automated External Defibrillator)

. . . . . . . .

## CONFIRM THAT THE PERSON CANNOT RESPOND

- Shout at and shake the person to make sure he or she isn't sleeping.
- Never shake an infant or young child.
   Instead, you can pinch the child to try to wake him or her up.

## CALL 9-1-1 OR HAVE SOMEONE ELSE CALL 9-1-1

#### TURN ON THE AED


- Make sure the person is in a dry location. If not, move them to a dry location.
- Most devices will give you step-by-step directions once turned on.

#### EXPOSE THE PERSON'S CHEST

• If the person's chest is wet, dry it.

#### PLACE THE AED PADS ON THE PERSON'S CHEST

- Place one pad on the right center of the person's chest above the nipple.
- Place the other pad slightly below the other nipple and to the left of the ribcage.

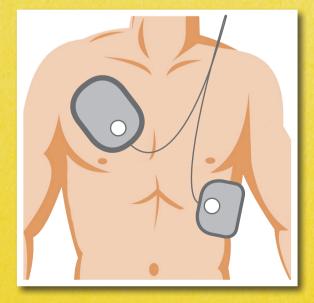

#### REMOVE METAL NECKLACES & UNDERWIRE BRAS

• The metal may conduct electricity & cause burns. You can cut the center of the bra and pull it away from the skin.

## CHECK THE PERSON FOR IMPLANTED MEDICAL DEVICES

- The outlines of these devices are visible under the skin on the chest or abdomen.
- · Also check for body piercings.
- Move the defibrillator pads at least 1 inch away from implated devices or piercings so the electric current can flow freely between the pads.

## FROM THE ELECTRODES ARE CONNECTED TO AED

- Make sure no one is touching the person, then press the AED's "analyze" button.
- Stay clear while the machine check's the person's heart rhythm.

### IF A SHOCK IS NEEDED, THE AED WILL TELL YOU WHEN TO DELIVER IT

 Stand clear of the person and make sure others are clear before you push the AED's "shock" button.

## START OR RESUME CPR UNTIL EMERGENCY MEDICAL HELP ARRIVES OR UNTIL PERSON BEGINS TO MOVE

 Stay with the peron until medical help arrives, and report all of the information you know about what has happened.







